Flowpression Boots Pro+

FLOWLIFE

Manual

FLOWLIFE


Active lifestyle. Made Possible.

Flowpression Boots Pro+

Premium Compression Boots by Flowlife

Flowpression Boots Pro+

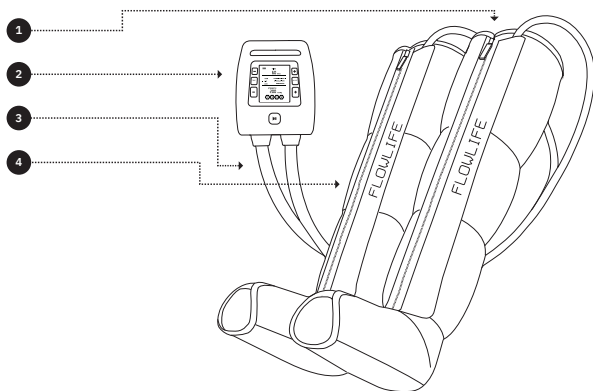1. Massage Boots
2. Control Unit
3. Air Hoses
4. Boot Zipper

6 — English
12 — Svenska
18 — Dansk

Flowpression Boots Pro+

Operating Instructions

EN

In this document Flowpression is referred as the product or Flowpression Boots Pro+.

The product must be used in accordance with the following safety instructions. This is to reduce the risk of fire, electric shock, property damage and personal injury. Read all safety instructions before using Flowpression.

SAFETY INSTRUCTIONS

SAFETY HAZARD - REDUCE RISK OF ELECTRIC SHOCK

- Never have the product plugged into an electrical outlet or other power source when cleaning the product.
- Never use the product near water, such as in a bath, shower, or similar.
- Never place the product anywhere where it risks being pulled into the water.
- Always store the product in a dry place, without contact with water or other liquids.
- Immediately unplug the product if it comes in contact with water.

DANGER - REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK AND PERSONAL INJURY

- Never leave the product unattended after connecting it to a power source.
- The product must (when used) never be covered with fabric or the like. This can lead to overheating and fire, electric shock and personal injury.
- With the conditions that the person in question has been informed about how the product is used in a safe way and about the safety risks stipulated in the safety instructions, children from the age of 13, people with various disabilities and people who do not have knowledge or experience of the product can use the product.
- Never use the product for any other purpose or in any other way than that specified in this safety instruction. The product should only be used for leg massage or compatible attachments for hips or arms. Never use product accessories from other than the original manufacturer.
- Never use the product in case of obvious or suspected damage to the product, the charging cord, the plug or the power source. If the product has come into contact with water, been dropped or damaged in any other way, you should always contact customer service for the possibility of inspection and repair. Never use the product until you have read the customer service information and recommendation.
- Never use the product cord as a handle to carry the product.
- Never store the product or cord on hot surfaces.

- Do not unzip the cuffs when using the product. Always turn off the power before unzipping the cuffs.
- Always be seated when using the product.
- Never insert anything into the product openings.
- Never use the product where aerosol products (spray products) are used or where oxygen is administered.
- Children must not use, clean or maintain the product without the supervision of an adult.
- Never try to repair a broken product on your own. Never open or unscrew the product. In the event of obvious or suspected damage to the product, customer service must always be contacted.
- If the product gets warm; turn it off and let it cool down before you continue using it.

CHARGE YOUR PRODUCT

- Flowpression should only be charged with the charger that comes with the product (12VDC/2A). Never leave product charging unattended. Never charge the product overnight. Connect and remove the charger carefully.

RISK GROUPS

If you belong to any of the following risk groups, ALWAYS consult your doctor before using Flowpression Boots Pro+. Only use the product with the approval of your doctor.

- Pregnancy, heart disease, medication *, diabetes with complications such as neuropathy or retinal damage, a pacemaker, recent surgery or injury, migraine or epilepsy, herniated disc, vertebral displacement, spondylosis or spondylolysis, recent arthroplasty or insertion of intrauterine inserts, metal pins or plates or other signs of concern about your physical health.

* *if you are taking medication of any kind, please read the package leaflet of your medicine before using Flowgression

RECYCLING AND SAFETY

- Flowpression contains a battery that must be safely returned to a recycling station
- Flowpressions battery is sorted as electronic waste.
- If the battery is exposed to extreme conditions or handled carelessly, it may pose a risk of fire or chemical burns. The battery should not be exposed to temperatures above 50 degrees Celsius. Never disassemble the battery. Never allow children to handle the battery.

INSTRUCTIONS FOR USE, MAINTENANCE AND STORAGE

To use the Flowpression Boots Pro+, follow the steps below:

- Find a comfortable position sitting, with your legs lying down or reclined.
- Press the power button on the back of the control unit.
- The Flowpression Boots Pro+ will start with sequence mode at 60 mmHg intensity on the default setting.
- Choose between strengths from 60-200mmHg, three different modes, and a zone treatment.
- You can choose between 3 program modes, Sequence, Circulation and Combination mode.
- Zone Mode: With the current Flowpression Boots Pro+, you can select zones specifically for the legs. However, by connecting the Hip or Arm attachment from Flowlife, you can access and use treatments for those areas as well.
- Press the -/+ intensity button to control the pressure between 60-200mmHg.
- Press the timer button to control the timer between 15-60 minutes.
- The product will stop its session when the time is up from when you set the timer.
- You can pause the treatment by pressing the PLAY/PAUS button.
- When you press the PLAY button, you only pause the unit. You need to turn off the device with the power button on the back and disconnect the hoses for complete deflation.
- Turn off the control unit by pressing the power button on the back

WARRANTY

If your product breaks within two years from the date of the purchase due to material or manufacturing defects, Flowlife will repair or replace necessary components or the product free of charge. Flowlife AB is not responsible for warranty and return matters for products purchased from another retailer.

The warranty is not valid if:

- The product is damaged by improper handling, misuse, accident or transport
- The product has been subjected to unauthorized repair
- The product has not been used according to the instructions from Flowlife
- The damage exceeds the cost of the product itself
- The product deteriorated after delivery due to abnormal storage conditions
- The dated proof of purchase has not been provided

Electrical specifications:

- Adapter
- AC Input: 100~240 Volts , AC 50/60Hz,
- DC Output: 12V 2A
- Rated power: 24W
- Compression Mode
- 3 Modes: Sequence / Circulation / Combination
- Intensity Levels / 60-200
- Pressure Range: 40-200mmhg
- Time: 15-60 minutes
- Product weight: 2.92 KG
- Package weight: 4.38 KG
- Sleeve Type: 4 - Chambers sleeve
- Foot / Lower Calf / Upper Calf / Thigh

- Battery: 2500mAh Battery Capacity
- 3-4 hours Use time
- 2.5 hours Charge time

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