Flowlight


FLOWLIFE

Manual

FLOWLIFE


Active lifestyle. Made Possible.

Flowlight

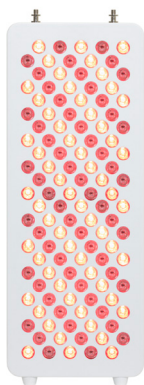Premium Red Light Therapy Panels by Flowlife

Flowlight

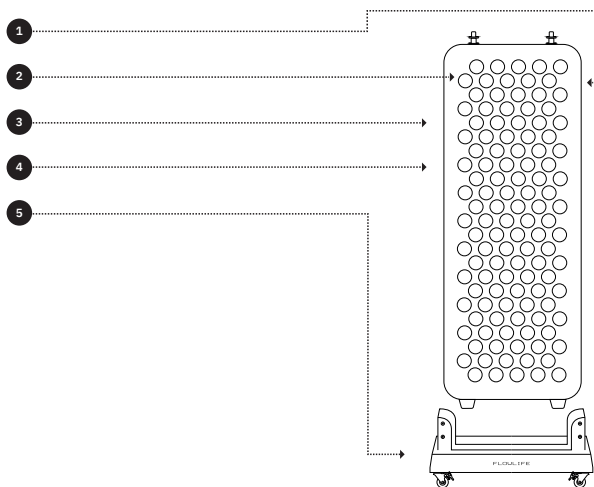1. Control Display
2. Infrared Chip Lamps
3. On/Of Switchch (On the backside)
4. Power Inlet (On the backside)
5. Floor Stand

6 — English
18 — Svenska

Flowlight

Operating Instructions

EN

Introduction

Welcome to your Flowlight Panel user manual. This guide will help you get started with your new red and near-infrared light therapy device and ensure you understand how to use it safely and effectively.

In this document, Flowlight is referred to as the product or Flowlight.

The following safety instructions must be used for the product. This reduces the risk of fire, electric shock, property damage, and personal injury. Read all safety instructions before using the Flowlight.

SAFETY INSTRUCTIONS

SAFETY HAZARD - REDUCE RISK OF ELECTRIC SHOCK

- Never have the product plugged into an electrical outlet or other power source when cleaning the product.
- Never use the product near water, such as in a bath, shower, or similar.
- Never place the product anywhere where it risks being pulled into the water.
- Always store the product in a dry place, without contact with water or other liquids.
- Immediately unplug the product if it comes in contact with water.

DANGER - REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK AND PERSONAL INJURY

- Never leave the product unattended after connecting it to a power source.
- The product must (when used) never be covered with fabric or the like. This can lead to overheating, fire, electric shock, and personal injury.
- Under the conditions that the person in question has been informed about how the product is used safely and about the safety risks stipulated in the safety instructions, children from the age of 15, people with various disabilities, and people who do not have knowledge or experience of the product can use it.
- Never use the product for any purpose other than that specified in this safety instruction. It should only be used for light therapy from a distance.
- Never use product accessories from other than the original manufacturer.
- Never use the product if there is apparent or suspected damage to it, the charging cord, the plug, or the power source.
- If the product has come into contact with water, been dropped, or damaged in any other way, you should always contact customer service for the possibility of inspection and repair.

- Only use the product once you have read the customer service information and recommendations.
- Never use the product cord as a handle to carry the product.
- Never store the product or cord on hot surfaces.
- Never insert anything into the product openings.
- Never use the product where aerosol products (spray products) are used or where oxygen is administered.
- Children may not use, clean, or maintain the product without the supervision of an adult.
- Never try to repair a broken product on your own, and never open or unscrew it. In the event of obvious or suspected damage to the product, customer service must always be contacted.
- Always consult your doctor before using the product if you need more clarification.

RISK GROUPS

If you belong to any of the following risk groups, ALWAYS consult your doctor before using Flowlight. Only use the product with the approval of your doctor.

- **Recent Burns:** Avoid using red light therapy on recent burns.
- **Heat Sensitivity:** Red light therapy generates heat, possibly posing risks in some conditions.
- **Light Sensitivity:** Individuals with extreme sensitivity to light or who have previously reacted to red light therapy should avoid its use.
- **Malignant Cancers:** There is a theoretical risk that red light therapy may increase cell proliferation in malignant cancers.
- **Hyperthyroidism:** Avoid red light therapy on the neck and upper chest if you have hyperthyroidism.
- **Pregnancy:** Avoid using red light therapy on the abdomen during pregnancy. Consult with a healthcare provider before use.
- **Eye Disease:** Consult with a healthcare provider if you have eye conditions before using red light therapy.
- **Fever or Infection:** Red light therapy may worsen existing hyperthermia caused by fever or infection.
- **Systemic Lupus Erythematosus (SLE):** Increased risk of reaction to red light therapy.
- **Severe Bleeding or Blood Loss:** Some research suggests red light therapy may slow the body's ability to stop bleeding.

- **Photosensitizing Medications:** Consult with a healthcare provider if you take medications that increase sensitivity to light.
- **Melasma/Hyperpigmentation:** Start slowly if you are concerned about hyperpigmentation.
- **Safety Glasses:** These are provided and should be worn when facing the lights.
- **Detox Symptoms:** This may cause detox-related symptoms initially. Start with shorter sessions and increase gradually.
- **Pre-existing Conditions/Medications:** Consult with a healthcare provider if you have pre-existing conditions or take photosensitizing drugs. Please consult your doctor if you have any other health issues or are unsure.

RECYCLING AND SAFETY

- Connect and remove the cord carefully.
- Flowlight should be returned to a recycling station as electronic waste.
- Avoid exposing Flowlight to extreme conditions or temperatures over 50°C to prevent the risk of fire.
- Never disassemble the Flowlight.
- Ensure the product is appropriately stable in the door hook or stand so it does not tip and fall.

WARRANTY

If your product breaks within two years from the date of the purchase due to material or manufacturing defects, Flowlife will repair or replace necessary components or the product free of charge. Flowlife AB is not responsible for warranty and return matters for products purchased from another retailer.

The warranty is not valid if:

- the product is damaged by improper handling, misuse, accident or transport
- the product has been subjected to unauthorized repair
- the product has not been used according to the instructions from Flowlife
- the damage exceeds the cost of the product itself
- the product deteriorated after delivery due to abnormal storage conditions
- the dated proof of purchase has not been provided

INSTRUCTIONS FOR USE AND MAINTENANCE

- Clean the light by unplugging it and wiping the surface with a cloth or paper towel. Avoid chemicals or bleach. Ensure the light is dry before using it again. Use indoors only.
- NIR LEDs: Near-infrared light is invisible to the naked eye. Use a camera to check if the LEDs are working.

INSTRUCTIONS FOR SETUP, USE AND STORAGE

Flowlight Panels Series Setup

In the package:

One panel. 1 Power cord. Manual. Remote Control. 2 pairs of eye protection glasses.

Start your Panel

- Connect the Flowlight to a power source and switch on the power switch on the back.

Programming

- You control everything through the LCD Screen on the side of your panel. The pause button allows you to start or stop the treatment. When you begin the session, the timer starts to count down.

Remote Control

- The remote control can also control all settings

Setting the timer

- Adjust the timer from 1 to 60 minutes based on your desired session length.

Adjusting light intensity

- Use the dimmer to set light intensity from 1% to 100%.

Pulse Setting

- Choose a pulse frequency between 1Hz to 10,000Hz.

During your session

- Relax and Enjoy: Sit back, relax, and let the Flowlight panel work.
- Monitor Settings: Adjust settings as needed for comfort and effectiveness.

Group: MASTER-SOLO-FOLLOWER

- When connecting multiple panels, you can connect up to four panels simultaneously.

WAVELENGTHS CONTROL

8 Wavelengths

- RED/NIR: Combines multiple wavelengths 610nm, 630nm, 660nm, 810nm, 830nm, 850nm, 940nm and 1060nm
- RED combines the wavelengths of 610nm, 630nm and 660nm.
- NIR combines the wavelengths of 810nm, 830nm, 850nm, 940nm and 1060nm.

7 Wavelengths

- RED/NIR: Combines multiple wavelengths 630nm, 660nm, 810nm, 830nm, 850nm, 940nm, and 1060nm.
- RED combines the wavelengths of 630nm and 660nm.
- NIR combines the wavelengths of 810nm, 830nm, 850nm, 940nm and 1060nm.

2 Wavelengths

- RED/NIR: Combines 660nm red light and 850nm near-infrared light.
- RED: Exclusively using the 660nm wavelength.
- NIR: Utilizing the 850nm wavelength.

ENDING YOUR SESSION

- Power Off: Turn off the device by pressing the power button. Unplug the device.
- Cool Down: Allow the panel to cool before handling or storing.
- Clean and Store: Wipe the panel with a soft, dry cloth. Store safely.

STANDARD USE

- First Time Users: Start with 10 minute sessions per area.
- General Use: We recommend 10-30 minute sessions 15-30 cm from your device. Allow at least 6 hours between sessions in the same area. Stay within 30 minutes in any one session. Use the device 4-12 times weekly, taking at least one day off.
- Eye Protection: Wear protective glasses when looking directly at the light, especially when NIR is active.

STORAGE

- Store in a cool, dry place. Unplug when not in use.

TECHNICAL SPECIFICATIONS

Flowlight 4300 Eight Waves

- Spectrum: 610, 630, 660, 810, 830, 850, 940, 1060nm
- Input voltage: 100-240V AC
- Frequency: 50/60Hz
- Lens: 30°
- Light density: $>130\text{mW/cm}^2$
- Dimension: 1660*505*65mm
- Digital Timer insert: default time is 30 minutes, Max 60mins
- Flicker: support 1~10000Hz;
- Dimmer: 1~100%
- Package included: 1pc light + 1pc EU plug + 1pc manual + 2pc glasses+ 1pc remote control
- Gross weight: 36kg/pc $\pm$ 1kg
- Rated power: 4320W
- Actual power: 1253W $\pm$ 10%
- 2.4G Remote Controller: Support

Flowlight 300 Two Waves

- Spectrum: 660nm, 850nm
- Input voltage: 100-240V AC
- Frequency: 50/60Hz
- Lens: 30°
- Light density: $>180\text{mW/cm}^2$
- Dimension: 290*210*65mm
- Default time: 30 minutes; max time is 60 minutes
- Pulse: support 1-10000Hz
- Dimmer: 1~100%
- Package included: 1pc light + 1pc wire with plug + 1pc manual + 2pc glasses+ 1pc Remote Controller
- Gross weight: 3.6kg/pc
- Rated power: 300W
- Actual power: About 120W
- Color: White

- Irradiance: About 168mw/cm² / 15 cm
- Power consumption: 0.12kWh

Flowlight 300 Seven Waves

- Spectrum: 630nm, 660nm, 810nm, 830nm, 850nm, 940nm, 1060nm
- Input voltage: 100 -240V AC
- Frequency: 50/60Hz
- Lens: 30°
- Light density: >180mW/cm²
- Dimension: 290*210*65mm
- Default time: 30 minutes; max time is 60 minutes
- Pulse: support 1-10000Hz
- Dimmer: 1~100%
- Package included: 1pc light + 1pc wire with plug + 1pc manual + 2pc glasses + 1pc Remote Controller
- Gross weight: 3.6kg/pc
- Rated power: 300W
- Actual power: About 120W
- Color: Black
- Irradiance: About 168mw/cm² / 15 cm
- Power consumption: 0.12kWh

Flowlight 600 Two Waves

- Spectrum: 660nm, 850nm
- Input voltage: 100 -240V AC
- Frequency: 50/60Hz
- Lens: 30°
- Light density: >130mW/cm²
- Dimension: 537*210*65mm
- Default time: 30 minutes; max time is 60 minutes
- Pulse: support 1-10000Hz;
- Dimmer: 1~100%
- Package included: 1pc light + 1pc wire with plug + 1pc manual + 2pc glasses + 1pc door hook + 1pc hanging cables + 1pc remote control
- Gross weight: 6.6kg/pc
- Rated power: 600W
- Actual power: About 250W
- Color: White
- Irradiance: About 176mw/cm² / 15 cm

- Power consumption: 0.25kWh

Flowlight 600 Seven Waves

- Spectrum: 630nm, 660nm, 810nm, 830nm, 850nm, 940nm, 1060nm
- Input voltage: 100 -240V AC
- Frequency: 50/60Hz
- Lens: 30°60°
- Light density: $>130\text{mW/cm}^2$
- Dimension: 537*210*65mm
- Default time: 30 minutes; max time is 60 minutes
- Pulse: support 1-10000Hz;
- Dimmer: 1~100%
- Package included: 1pc light + 1pc wire with plug + 1pc manual + 2pc glasses + 1pc door hook + 1pc hanging cables + 1pc remote control
- Gross weight: 6.6kg/pc
- Rated power: 600W
- Actual power: About 250W
- Color: Black
- Irradiance: About 176mw/cm^2 / 15 cm
- Power consumption: 0.25kWh

Flowlight 1500 Two Waves

- Spectrum: 660nm, 850nm
- Input voltage: 100 -240V AC
- Frequency: 50/60Hz
- Lens: 30°
- Light density: $>130\text{mW/cm}^2$
- Dimension: 920*300*65mm
- Default time: 30 minutes; max time is 60 minutes
- Pulse: support 1-10000Hz;
- Dimmer: 1~100%
- Package included: 1pc light + 1pc wire with plug + 1pc manual + 2pc glasses + 1pc door hook + 1pc hanging cables + 1pc remote control
- Gross weight: 13.5kg/pc
- Rated power: 1500W
- Actual power: About 540W
- Color: White
- Irradiance: About 190mw/cm^2 / 15 cm
- Power consumption: 0.55kWh

Flowlight 1500 Seven Waves

- Spectrum: 630nm, 660nm, 810nm, 830nm, 850nm, 940nm, 1060nm
- Input voltage: 100 - 240V AC
- Frequency: 50/60Hz
- Lens: 30°60°
- Light density: $>130\text{mW/cm}^2$
- Dimension: 920*300*65mm
- Default time: 30 minutes; max time is 60 minutes
- Pulse: support 1-10000Hz;
- Dimmer: 1~100%
- Package included: 1pc light + 1pc wire with plug + 1pc manual + 2pc glasses + 1pc door hook + 1pc hanging cables + 1pc remote control
- Gross weight: 13.5kg/pc
- Rated power: 1500W
- Actual power: About 540W
- Color: Black
- Irradiance: About 190mw/cm^2 / 15 cm
- Power consumption: 0.55kWhWERWER

STAND ASSEMBLY INSTRUCTIONS

****Warning:**** Adult assembly is required. Keep children away during assembly. This product contains small parts that may present a choking hazard to children under three years old.

Tilt risk. This panel is significant and is at the risk of tilting if not used and handled correctly. Please read the Manual before use.


- Don't push or lean against the panel.
- Turn off the power switch and unplug the power cord before you move the panel.
- Unlock the wheels if they are locked before you move it.
- When moving it, use the wheels and two hands on each side of the panels' sides to make sure you can hold on to it if something happens and would make a risk of tilting the panel.
- Make sure no children are around when you move the panel.
- Make sure no children are close to the panel without adult supervision.
- Make sure no activities around the panel can pose a risk of forcing it to tilt.
- Keep a distance from the panel when it's not in use.

****Using the Floor Stand****



To set up your panel with the floor stand, follow these steps:

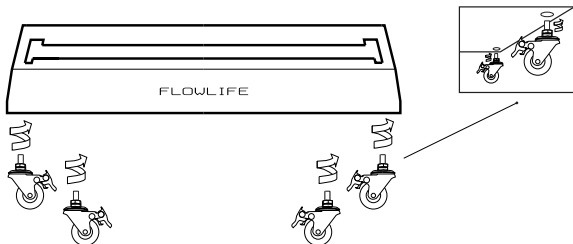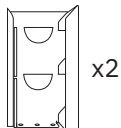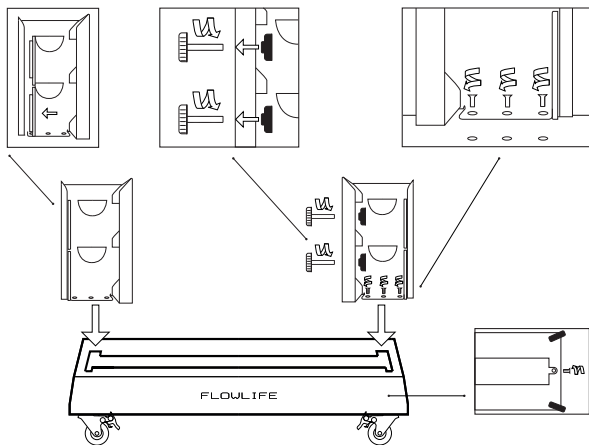1. Press the black plate on one of the wheels to lock its function.
2. Flip the stand base upside down.
3. Screw the wheels into the four corners of the base.
4. Insert the four soft foam lines into each bracket.
5. Unscrew the six screws from each side of the stand's pathway.
6. Secure the brackets by screwing the screws back into the stand through the mounting brackets, using the same holes from which you removed the screws. Ensure the soft foam lines are correctly aligned and do not twist while fastening the mounting brackets.
7. Secure and fasten the pathway by screwing two screws on each side of the brackets from the bottom.
8. Screw the wheel screws from the back of the brackets, leaving approximately half of the screw length exposed.
9. Attach the black protection pads on top of the screws, then screw the wheel screws back in to create space for the panel.
10. Attach the panel to the stand, ensuring that the wheels are locked and there is enough space in the pathway between the brackets.
11. Tighten the black wheel screws until they feel firm.

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