

KEEP YOUR PLAN IN PLACE

Reliable, Predictable, Pliable and Fully Resorbable.



 **conFORM**

resorbable collagen membrane

- Fully resorbable with a predictable resorption time of 3-4 months¹
- Easy to handle, conforms to the surgical site without re-shaping and does not adhere to itself
- Guides healing and regeneration of bone and surrounding tissue
- Unique fiber orientation provides high tensile strength for *in situ* stability of membrane
- Cell occlusive, designed to retard and/or prevent epithelium down growth

¹ Reference on back



conFORM RESORBABLE COLLAGEN MEMBRANE

conFORM is a resorbable collagen membrane matrix made from Type I Bovine Collagen. Its pliable characteristics and predictable resorption time of 3-4 months make it an ideal membrane for guiding the healing of bone and surrounding tissue. When hydrated, conFORM will naturally adapt to surgical defects and maintain a high tensile strength, which makes it the ideal membrane for suturing.

With its in vivo stability, permeability and mechanical strength, conFORM is a reliable choice for various clinical applications. When used in conjunction with a bone grafting particulate, conFORM can help mitigate the migration of the particulate from the surgical site. Its resorbability eliminates the need for surgical removal. conFORM is packaged sterile, non-pyrogenic, and for single use.

Features

- White, non-friable membrane matrix engineered from highly purified Type I Collagen
- Fully resorbable with a predictable resorption time of 3-4 months ¹
- Guides healing and regeneration of bone and surrounding tissue
- Unique fiber orientation provides high tensile strength for *in situ* stability of membrane
- Cell occlusive, designed to retard and/or prevent epithelium down growth

Advantages


- Conforms to surgical site without re-shaping
- Easy to handle, does not adhere to itself
- Sufficient resorption time for sustained function in bone repair
- Can easily be trimmed during surgery for a final fit
- May be placed with either side against the tissue or bone

Typical Uses

- Ridge augmentation
- Bone regeneration around implants
- Sinus lift
- Periodontal defects
- Apical root resection
- Extraction sites



resorbable collagen membrane

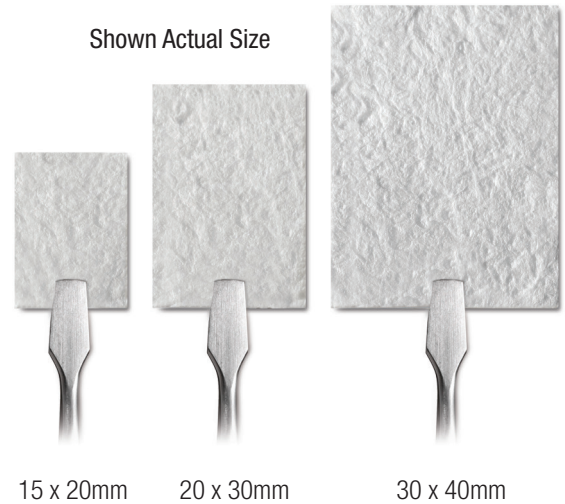

conFORM			Regular	Member	Premier	Premier+
5091520	15 x 20mm	each	\$146.99	\$146.99	\$124.94	\$124.94
5092030	20 x 30mm	each	\$181.29	\$181.29	\$154.10	\$154.10
5093040	30 x 40mm	each	\$291.99	\$291.99	\$248.19	\$248.19

Copyright © 2023 ACE SOUTHERN

All prices are subject to change without notice. Prices shown are domestic prices in U.S. dollars; does not include shipping or any duties. For discount details visit acesouthern.com/ASAP. And for the most up-to-date pricing, please contact ACE SOUTHERN or visit our website.

Footnote References:

- ¹ D.Yuen, J. Ulreich, G. Zuclich, H-B. Lin and S-T Li "Prediction of In Vivo Stability of a Resorbable, Reconstituted Type I Collagen Membrane by In Vitro Methods." 2000 Society for Biomaterials - Sixth World Congress.