

# SCREW IT | TACK IT | SECURE IT

Unmatched Versatility



complete fixation system

- Everything you need for both bone fixation and membrane stabilization
- Incorporates all the components necessary to pick-up and drive truTACKS and truSCREWS
- True centered locking connection allows for axial alignment with each pick-up
- Removal sleeve for disengaging the screw from the driver, without damaging screw head
- truSCREWS and truTACKs are developed for easy insertion with maximum fixation



## complete fixation system

The truFIX System from ACE provides you with everything you need for both bone fixation and membrane stabilization. This system incorporates all the necessary components to pick-up and drive the ACE truSCREW and truTACK. The truFIX system eliminates the need for multiple systems and unnecessary components, making it user-friendly for your practice.

The truSCREW, with its aggressive cutting flutes, is the ideal bone screw for the fixation of small bone within the oral and maxillofacial environment. These cutting flutes (in most instances) eliminate the need for any pre-drilling. The design of the screw ensures an effortless insertion into all types of bone.

The truTACK makes the stabilization of membranes quick and trouble-free. ACE's unique tack incorporates a hexagon on its head and threads on its shaft, allowing for easy removal. The truTACK is placed like a tack and removed like a screw, a feature that you will not find in any other system.

The truTENT is a refinement of the truSCREW. Its raised collar and wider head is designed to support a titanium mesh during augmentation procedures.

- True centered locking connection
- True axial alignment with pickup each time
- truSCREW self-drilling screw developed for easy insertion with maximum fixation
- truSCREW removal sleeve for disengaging the screw from the driver without damaging screw head
- truTACK hexagon driver head and barbed tip to pierce without drilling for simple tack insertion and easy screw-like removal



High strength surgical grade titanium alloy (6Al-4V ELI)


Hex head driver

Smooth beveled seating neck


Thread for easier removal

Anchor collar to grab and secure tack in place

Barbed tip to pierce without pre-drilling

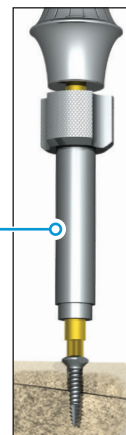


The driver's true-centered locking connection holds the truSCREW securely.

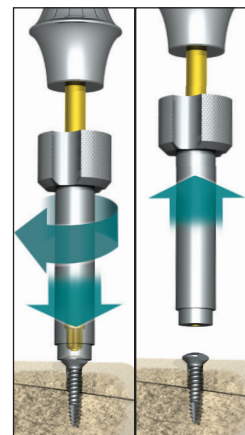


The removal sleeve disengages the driver from the screw without compromising the screw or the stability of the graft.



**D** truSCREW Driver Tip **C** and truSCREW Driver Removal Sleeve **D** are shown together as they come in the kit. Replacements are sold individually.

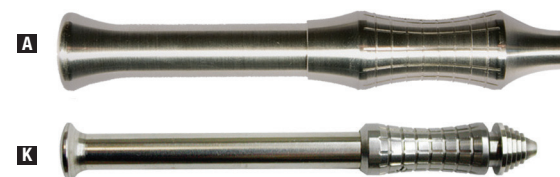

Optional Small Driver Handle shown next to the truFIX Driver Handle to show size difference

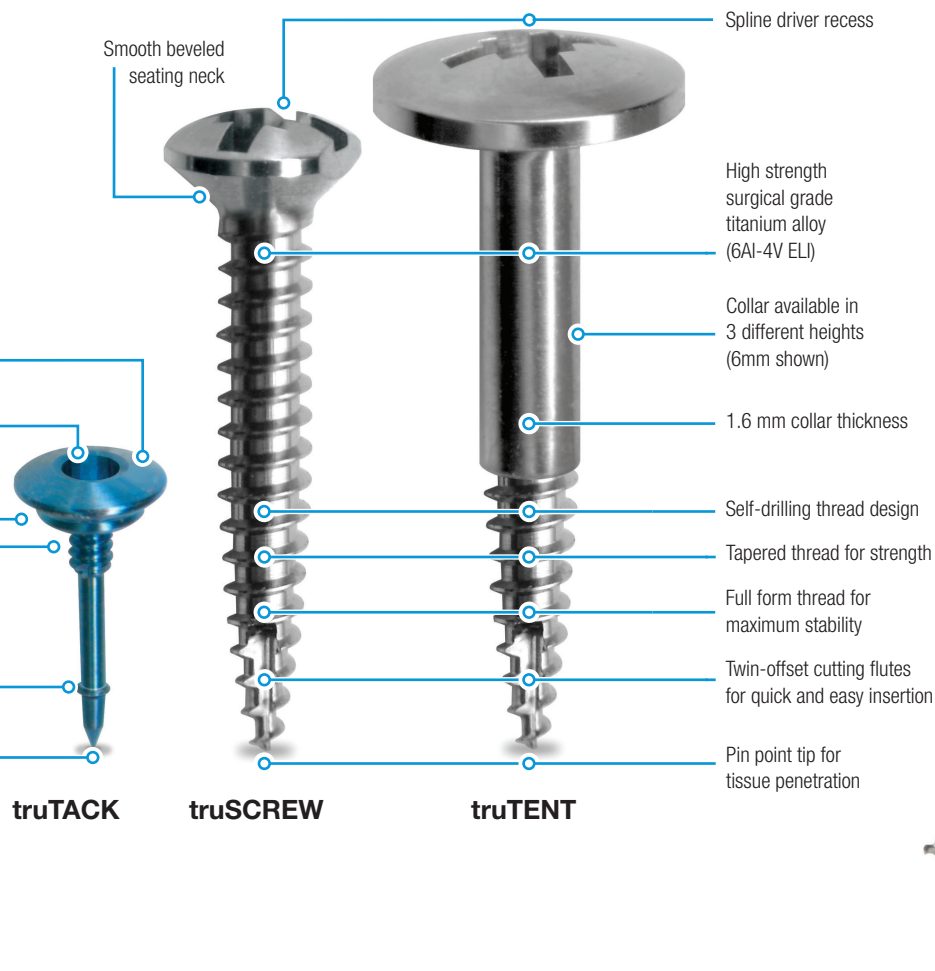

## truFIX KIT AND COMPONENTS

			Regular	Member	Premier	Premier+
45418015	<b>truFIX Complete Fixation System</b>	1 kit	<b>\$1,326.79</b>	<b>\$1,326.79</b>	<b>\$1,194.11</b>	<b>\$1,194.11</b>

### System Includes:

45419001	<b>A</b> truFIX Driver Handle	each	<b>\$372.29</b>	<b>\$372.29</b>	<b>\$335.06</b>	<b>\$335.06</b>
45417001	<b>B</b> truTACK Driver Tip (comes with Tip Cover)	each	<b>\$113.29</b>	<b>\$113.29</b>	<b>\$101.96</b>	<b>\$101.96</b>
45415001	<b>C</b> truSCREW Driver Tip	each	<b>\$88.29</b>	<b>\$88.29</b>	<b>\$79.46</b>	<b>\$79.46</b>
45417901	<b>D</b> truSCREW Driver Removal Sleeve	each	<b>\$91.79</b>	<b>\$91.79</b>	<b>\$82.61</b>	<b>\$82.61</b>
45415201	<b>E</b> truSCREW Driver Tip, Contra-Angle	each	<b>\$100.29</b>	<b>\$100.29</b>	<b>\$90.26</b>	<b>\$90.26</b>
45440203	<b>F</b> 2-Step Countersink Bur, Ø 0.8 and 1.6mm, Contra-Angle	each	<b>\$72.99</b>	<b>\$72.99</b>	<b>\$65.69</b>	<b>\$65.69</b>
45440202	<b>G</b> Ø 1.1mm Pilot Twist Drill, Contra-Angle, 29.8mm Length	each	<b>\$84.79</b>	<b>\$84.79</b>	<b>\$76.31</b>	<b>\$76.31</b>
502700045	<b>H</b> Ø 0.45mm Pilot Bur, Contra-Angle, 27mm Length	each	<b>\$55.79</b>	<b>\$55.79</b>	<b>\$50.21</b>	<b>\$50.21</b>
20157702	<b>I</b> Hand Driver, 0.88mm External Hex	each	<b>\$87.49</b>	<b>\$87.49</b>	<b>\$78.74</b>	<b>\$78.74</b>
45450201	<b>J</b> truSCREW Packaging Removal Tool	each	<b>\$55.79</b>	<b>\$55.79</b>	<b>\$50.21</b>	<b>\$50.21</b>
45418501	truFIX Tray	each	<b>\$435.49</b>	<b>\$435.49</b>	<b>\$391.94</b>	<b>\$391.94</b>

### truFIX Small Driver Handle (sold separately)

This smaller, thinner, and lighter handle is often preferred by surgeons with smaller hands.

45419501	<b>K</b> Small Driver Handle, 89mm (3.5") long	each	<b>\$148.49</b>	<b>\$148.49</b>	<b>\$133.64</b>	<b>\$133.64</b>
----------	--	------	-----------------	-----------------	-----------------	-----------------



## STERILE TACKS AND SCREWS

truTACKS		Head Ø	Thread Ø x Length		Regular	Member	Premier	Premier+
9600313	truTACK (gray)	2.5mm	0.7mm x 3.0mm	10/pkg	<b>\$167.99</b>	<b>\$167.99</b>	<b>\$142.79</b>	<b>\$142.79</b>
9600314	truTACK (blue)	2.5mm	0.7mm x 5.0mm	10/pkg	<b>\$167.99</b>	<b>\$167.99</b>	<b>\$142.79</b>	<b>\$142.79</b>

### truSCREW 1.2 Thread Diameter

45427202	truSCREW	3mm HD	1.2mm x 4.5mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45427203	truSCREW	3mm HD	1.2mm x 6.0mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45427204	truSCREW	3mm HD	1.2mm x 7.5mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45427205	truSCREW	3mm HD	1.2mm x 9.0mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45427206	truSCREW	3mm HD	1.2mm x 10.5mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>

### truSCREW 1.5 Thread Diameter

45427502	truSCREW	3mm HD	1.5mm x 6.0mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45427503	truSCREW	3mm HD	1.5mm x 7.5mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45427504	truSCREW	3mm HD	1.5mm x 9.0mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45427505	truSCREW	3mm HD	1.5mm x 10.5mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45427506	truSCREW	3mm HD	1.5mm x 12.0mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45427507	truSCREW	3mm HD	1.5mm x 13.5mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45427508	truSCREW	3mm HD	1.5mm x 15.0mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>

### truSCREW 2.0 Thread Diameter

45428002	truSCREW	3mm HD	2.0mm x 6.0mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45428003	truSCREW	3mm HD	2.0mm x 7.5mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45428004	truSCREW	3mm HD	2.0mm x 9.0mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45428005	truSCREW	3mm HD	2.0mm x 10.5mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45428006	truSCREW	3mm HD	2.0mm x 12.0mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45428007	truSCREW	3mm HD	2.0mm x 13.5mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>
45428008	truSCREW	3mm HD	2.0mm x 15.0mm	5/pkg	<b>\$163.29</b>	<b>\$163.29</b>	<b>\$138.80</b>	<b>\$138.80</b>

### truTENT TENTING SCREW 1.5 Thread Diameter

454391001	truTENT 4mm collar	5mm HD	1.5mm x 10.0mm	each	<b>\$59.79</b>	<b>\$59.79</b>	<b>\$50.82</b>	<b>\$50.82</b>
454391002	truTENT 6mm collar	5mm HD	1.5mm x 12.0mm	each	<b>\$59.79</b>	<b>\$59.79</b>	<b>\$50.82</b>	<b>\$50.82</b>
454391003	truTENT 8mm collar	5mm HD	1.5mm x 14.0mm	each	<b>\$59.79</b>	<b>\$59.79</b>	<b>\$50.82</b>	<b>\$50.82</b>

Copyright © 2024 ACE SOUTHERN

All prices are subject to change without notice. Prices shown are domestic prices in U.S. dollars; does not include shipping or any duties. For discount details visit [acesouthern.com](https://acesouthern.com)/ASAP. And for the most up-to-date pricing, please contact ACE SOUTHERN or visit our website.