



ACE SOUTHERN

A HENRY SCHEIN® BRAND

# SCREW IT | TACK IT | SECURE IT

Unmatched Versatility



complete fixation system

- Everything you need for both bone fixation and membrane stabilization
- Incorporates all the components necessary to pick-up and drive truTACKS and truSCREWS
- True centered locking connection allows for axial alignment with each pick-up
- Removal sleeve for disengaging the screw from the driver, without damaging screw head
- truSCREWS and truTACKs are developed for easy insertion with maximum fixation



## complete fixation system


The truFIX System from ACE provides you with everything you need for both bone fixation and membrane stabilization. This system incorporates all the necessary components to pick-up and drive the ACE truSCREW and truTACK. The truFIX system eliminates the need for multiple systems and unnecessary components, making it user-friendly for your practice.

The truSCREW, with its aggressive cutting flutes, is the ideal bone screw for the fixation of small bone within the oral and maxillofacial environment. These cutting flutes (in most instances) eliminate the need for any pre-drilling. The design of the screw ensures an effortless insertion into all types of bone.

The truTACK makes the stabilization of membranes quick and trouble-free. ACE's unique tack incorporates a hexagon on its head and threads on its shaft, allowing for easy removal. The truTACK is placed like a tack and removed like a screw, a feature that you will not find in any other system.

The truTENT is a refinement of the truSCREW. Its raised collar and wider head is designed to support a titanium mesh during augmentation procedures.

- True centered locking connection
- True axial alignment with pickup each time
- truSCREW self-drilling screw developed for easy insertion with maximum fixation
- truSCREW removal sleeve for disengaging the screw from the driver without damaging screw head
- truTACK hexagon driver head and barbed tip to pierce without drilling for simple tack insertion and easy screw-like removal

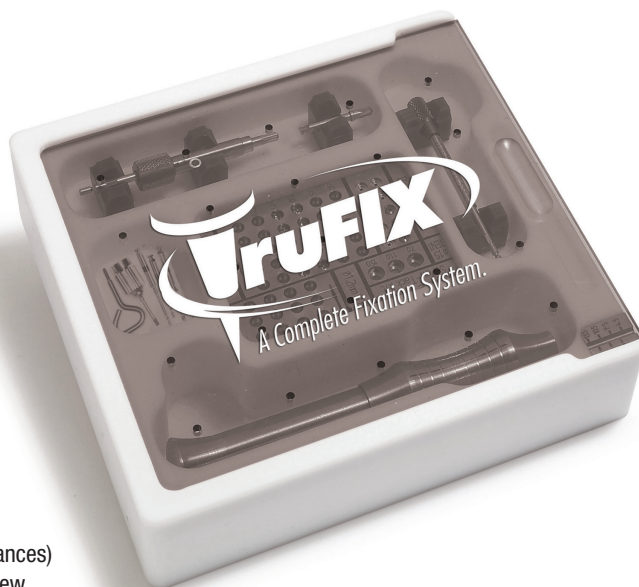

High strength surgical grade titanium alloy (6Al-4V ELI)


Hex head driver

Smooth beveled seating neck


Thread for easier removal

Anchor collar to grab and secure tack in place

Barbed tip to pierce without pre-drilling

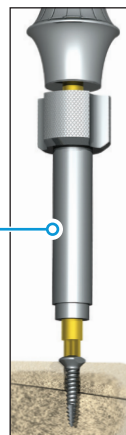

The driver's true-centered locking connection holds the truSCREW securely.

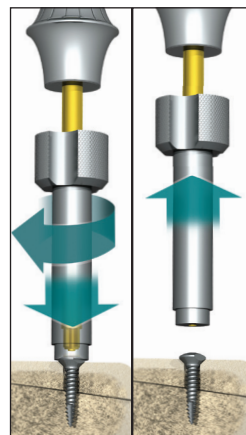


The removal sleeve disengages the driver from the screw without compromising the screw or the stability of the graft.



The blue protective Tip Cover is sold together with the truTACK Driver Tip **B**.

**D** truSCREW Driver Tip **C** and truSCREW Driver Removal Sleeve **D** are shown together as they come in the kit. Replacements are sold individually.

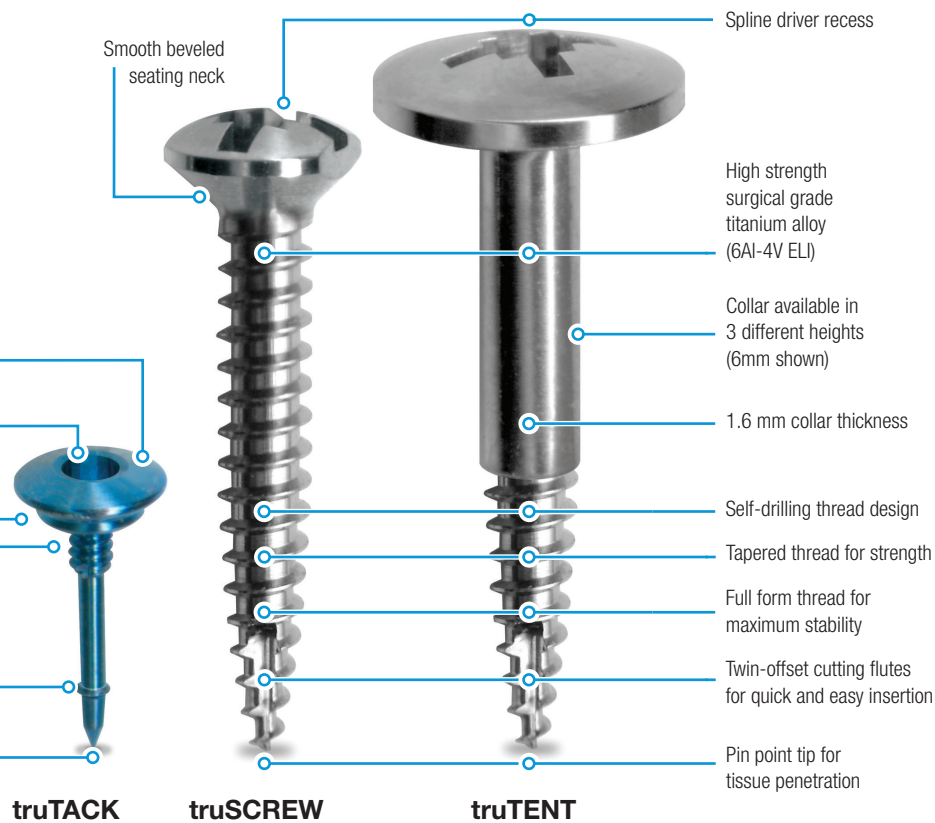

## truFIX KIT AND COMPONENTS

1101796      **truFIX Complete Fixation System**      each

### System Includes:

1464523	<b>A</b>	truFIX Driver Handle	each
1464535	<b>B</b>	truTACK Driver Tip (comes with Tip Cover)	each
1464534	<b>C</b>	truSCREW Driver Tip	each
1476299	<b>D</b>	truSCREW Driver Removal Sleeve	each
1464533	<b>E</b>	truSCREW Driver Tip, Contra-Angle	each
1476166	<b>F</b>	2-Step Countersink Bur, Ø 0.8 and 1.6mm, Contra-Angle	each
1467036	<b>G</b>	Ø 1.1mm Pilot Twist Drill, Contra-Angle, 29.8mm Length	each
1466789	<b>H</b>	Ø 0.45mm Pilot Bur, Contra-Angle, 27mm Length	each
1469232	<b>I</b>	Hand Driver, 0.88mm External Hex	each
1476167	<b>J</b>	truSCREW Packaging Removal Tool	each
truFIX Tray (not sold individually)			



## STERILE TACKS AND SCREWS

truTACK		Thread Ø x Length	Head Ø	
1107029	gray	0.7mm x 3.0mm	2.5mm	10/pkg
1102392	blue	0.7mm x 5.0mm	2.5mm	10/pkg

### truSCREW 1.2 Thread Diameter

1108879		1.2mm x 4.5mm	3mm	5/pkg
1102586		1.2mm x 6.0mm	3mm	5/pkg
1109615		1.2mm x 7.5mm	3mm	5/pkg
1104330		1.2mm x 9.0mm	3mm	5/pkg
1105857		1.2mm x 10.5mm	3mm	5/pkg

### truSCREW 1.5 Thread Diameter

1107821		1.5mm x 6.0mm	3mm	5/pkg
1105111		1.5mm x 7.5mm	3mm	5/pkg
1105979		1.5mm x 9.0mm	3mm	5/pkg
1103433		1.5mm x 10.5mm	3mm	5/pkg
1104087		1.5mm x 12.0mm	3mm	5/pkg
1108365		1.5mm x 13.5mm	3mm	5/pkg
1102598		1.5mm x 15.0mm	3mm	5/pkg

### truSCREW 2.0 Thread Diameter

1467077		2.0mm x 6.0mm	3mm	5/pkg
1467075		2.0mm x 7.5mm	3mm	5/pkg
1467076		2.0mm x 9.0mm	3mm	5/pkg
1467073		2.0mm x 10.5mm	3mm	5/pkg
1467074		2.0mm x 12.0mm	3mm	5/pkg
1467071		2.0mm x 13.5mm	3mm	5/pkg
1467072		2.0mm x 15.0mm	3mm	5/pkg

### truTENT Tenting Screw 1.5 Thread Diameter

1464486	4mm collar	1.5mm x 10.0mm	5mm	each
1464447	6mm collar	1.5mm x 12.0mm	5mm	each
1467078	8mm collar	1.5mm x 14.0mm	5mm	each

Copyright © Henry Schein, Inc.

Henry Schein is a registered trademark of HS TM, LLC.

For the most up-to-date pricing, please contact ACE SOUTHERN or visit our website at [www.acesouthern.com](http://www.acesouthern.com)